



Centre Anti-Poison pour le Québec: (800) 463-5060

Tél. (Qc): (418) 660-8666 / 800-890-8666

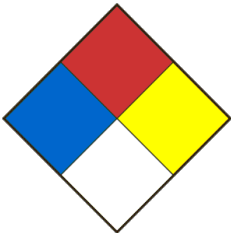
Fax. (Qc): (418) 660-8998

SAFETY DATA SHEET

SECTION 01 - PRODUCT AND COMPANY IDENTIFICATION

Product Identifier POTASSIUM CHROMATE 10%W/V		Product Use Laboratory use	
Chemical formula K_2CrO_4		Product code PS-0887	Molar weight 194,2
Chemical name / Commercial name / Synonymous POTASSIUM CHROMATE, POTASSIUM CHROMATE (VI), NEUTRAL POTASSIUM CHROMATE, BIPOTASSIUM CHROMATE, DIPOTASSIUM CHROMATE, DIPOTASSIUM MONOCHROMATE			
Supplier's name Laboratoire MAT		Address-Street 610, Adanac Street	
City Québec		Province Québec	
Postal code G1C 7B7	Internet www.labmat.com	Phone number 418-660-8666 / 800-890-8666	
Emergency phone	CANUTEC: 613-996-6666 CENTRE ANTI-POISON DU QUÉBEC 800-463-5060		
Date SDS 4/14/2023	SDS Prepared by Laboratoire MAT	E-Mail labmat@labmat.com	

SECTION 02 - HAZARDS IDENTIFICATION

Classification WHIMS / GHS	Skin corrosion/irritation - Skin irritation category 2 Serious eye damage/eye irritation - Serious eye damage category 1 Respiratory or skin sensitization - Respiratory sensitizer category 1 Acute toxicity - Oral category 4 Germ cell mutagenicity category 1 Carcinogenicity category 1 Respiratory or skin sensitization - Skin sensitize category 1	
Signal Word	DANGER	
Hazards statements (H)	H340 May cause genetic defects. H350 May cause cancer. H302 Harmful if swallowed. H315 Causes skin irritation. H318 Causes serious eye damage. H334 May cause allergy or asthma symptoms or breathing difficulties if inhaled. H317 May cause an allergic skin reaction.	
Precautionary statements (P)	P201 Obtain special instructions before use. P202 Do not handle until all safety precautions have been read and understood. P261 Avoid breathing dust / fume / gas / mist / vapors / spray. P280 Wear protective gloves/protective clothing/eye protection/face protection. P302 + P352 IF ON SKIN: Wash with plenty of soap and water. P308 + P313 IF exposed or concerned: Get medical advice/attention. P321 Specific treatment (see section 4 of the SDS and on this label). P362 + P364 Take off contaminated clothing and wash it before reuse. P405 Store locked up. P501 Dispose of contents/container in accordance with local / regional / national / international regulations or contact a specialist waste disposal company. P264 Wash the areas of the body that have been in contact with the product after handling. P270 Do no eat, drink or smoke when using this product. P284 Wear respiratory protection. P301 + P312 IF SWALLOWED: Call a POISON CENTER or doctor/physician if you feel unwell. P304 + P340 IF INHALED: Remove victim to fresh air and keep at rest in a position comfortable for breathing. P305 + P351 + P338 IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. P310 Immediately call a POISON CENTER or doctor/physician. P330 Rinse mouth. P332 + P313 If skin irritation occurs: Get medical advice/attention. P342 + P311 If experiencing respiratory symptoms: Call a POISON CENTER or doctor/physician. P272 Contaminated work clothing should not be allowed out of the workplace. P333 + P313 If skin irritation or rash occurs: Get medical advice/attention.	
PICTOGRAMS		
Other dangers	NFPA (Risk: 0=No risk; 1=Slight; 2=Moderate; 3=Signifiant; 4=Extreme)	
	Health 3 Fire 0 Reactivity 0 Special danger	

SECTION 03 - COMPOSITION/INFORMATION ON INGREDIENTS

Ingrédients (Dénomination chimique / synonymes)	Numéro CAS et tout identificateur unique	Concentration (%)
Chromate de potassium	7789-00-6	9

SECTION 04 - FIRST AID MEASURES

Eye contact	Wash eyes with large amounts of water for at least 15 minutes while holding eyelids apart to rinse eyes. If irritation persists, seek medical attention.
Skin contact	Wash skin with plenty of water for at least 15 minutes. Remove soiled clothing. If irritation persists, seek medical attention.
Inhalation	Move the unwell person to the fresh air. If breathing is difficult, give oxygen. Consult a physician.
Ingestion	If the person is conscious, give water to drink. Never give anything by mouth to an unconscious person. Consult a physician.
Most important symptoms and effects (acute and delayed)	The product is toxic, corrosive and irritant. Main symptoms of high exposure: Skin, eye and respiratory system irritation. Skin sensitizer. Vomiting. Diarrhea. Hemorrhage. Liver damage. Damage to the digestive system. Runny nose. Sneezing. Ref. section 11.
Immediate medical attention and special treatment, if necessary	In case of medical consultation, keep this sheet available.
General advice	Show this safety data sheet to the doctor in attendance.

SECTION 05 - FIREFIGHTING MEASURES

Flammability	No
Ignition conditions	Not flammable or combustible. May cause fire on contact with organic material. May ignite if product dries on clothing, wood or other combustible materials.
Suitable extinguishing media	Use water to extinguish the fire
Unsuitable extinguishing media	CO2 or halon can provide limited control.
Hazardous combustion products	Hazardous combustion products formed under fire conditions: - Potassium oxides, Chromium oxides.
Special fire and explosion hazards	Oxidizer. In the event of a fire, the water in the solution will evaporate, and the oxidizing base product will fuel the fire. May react violently with incompatible products (Ref Section 10).
Special protective equipment and precautions for firefighters	Wear self contained breathing apparatus for fire fighting if necessary.

SECTION 06 - ACCIDENTAL RELEASE MEASURES

Methods and materials for containment and cleaning up / Personnel precautions, protective equipment	Evacuate personnel to safe areas. Absorb the product with sand or vermiculite. Do not use cellulose-based absorbent. Dilute residues with water, clean and rinse. Ensure a good ventilation of the premises. When handling, wear appropriate safety equipment. Use a respirator as needed. Dispose of residues in a container provided for the disposal of hazardous materials. Discharge into the environment must be avoided.
--	---

SECTION 07 - HANDLING AND STORAGE

Conditions for safe storage	Store in a cool, dry place. Keep container tightly closed and store away from heat, moisture, and incompatible products.
Methods of handling	Avoid contact with the skin, eyes and clothes. Avoid contact with combustible materials. Provide appropriate exhaust ventilation at places where dust or vapor is formed. Avoid ingestion and inhalation. Wear personal protective equipment when handling. Always ensure good ventilation. Transport according to TDG (ref Section 14)

SECTION 08 - EXPOSURE CONTROLS/PERSONAL PROTECTION

Workplace control parameters

Components	CAS-No.	Value	Control parameters	Basis
Potassium chromate	7789-00-6			Canada. Alberta, Occupational Health and Safety Code (table 2: OEL)
		TWA	0.050000 mg/m3	Canada. Alberta, Occupational Health and Safety Code (table 2: OEL)
Remarks	Confirmed Human Carcinogen (means that the agent is carcinogenic to humans)			
		TWAEV	0.050000 mg/m3	Québec. Regulation respecting occupational health and safety, Schedule 1, Part 1: Permissible exposure values for airborne contaminants
	A substance which may not be recirculated in accordance with section 108 A substance to which exposure must be reduced to a minimum in accordance with section 42 Sensitizer Carcinogenic effect detected in humans			
		TWA	0.025000 mg/m3	Canada. British Columbia OEL
	IARC '1' applies to substances categorized as carcinogenic to humans, and used when there is sufficient evidence of carcinogenicity in humans. ACGIH 'A1' applies to those substances confirmed as human carcinogens based on the weight of evidence from epidemiological studies			
		C	0.100000 mg/m3	Canada. British Columbia OEL
	IARC '1' applies to substances categorized as carcinogenic to humans, and used when there is sufficient evidence of carcinogenicity in humans. ACGIH 'A1' applies to those substances confirmed as human carcinogens based on the weight of evidence from epidemiological studies			
				Canada. Alberta, Occupational Health and Safety Code (table 2: OEL)
		TWA	0.05 mg/m3	Canada. Alberta, Occupational Health and Safety Code (table 2: OEL)
	Confirmed Human Carcinogen (means that the agent is carcinogenic to humans)			
		TWAEV	0.05 mg/m3	Québec. Regulation respecting occupational health and safety, Schedule 1, Part 1: Permissible exposure values for airborne contaminants
	A substance which may not be recirculated in accordance with section 108 A substance to which exposure must be reduced to a minimum in accordance with section 42 Sensitizer Carcinogenic effect detected in humans			
		TWA	0.025 mg/m3	Canada. British Columbia OEL
	IARC '1' applies to substances categorized as carcinogenic to humans, and used when there is sufficient evidence of carcinogenicity in humans. ACGIH 'A1' applies to those substances confirmed as human carcinogens based on the weight of evidence from epidemiological studies			
		C	0.1 mg/m3	Canada. British Columbia OEL
	IARC '1' applies to substances categorized as carcinogenic to humans, and used when there is sufficient evidence of carcinogenicity in humans. ACGIH 'A1' applies to those substances confirmed as human carcinogens based on the weight of evidence from epidemiological studies			

Data source	Sigma-Aldrich.
Ventilation	Fan.
Respiratory	If work under the hood is not possible, or if the permissible levels are exceeded, use a mechanical filter / cartridge against NIOSH vapors or a respirator with air supply.
Gloves	Handle with gloves.
Eyes	Safety goggles with safety shutters.
Shoes	Safety shoes.
Clothing	Labcoat. The type of protective equipment must be selected according to the concentration and amount of the dangerous substance at the specific workplace.
Engineering control	Have safety showers and eyewash stations in the workplace in case of an emergency and a ventilation system to maintain the level of concentrations in the air below the exposure limit values.

SECTION 09 - PHYSICAL AND CHEMICAL PROPERTIES

Physical state	Liquid.
Appearance	Liquide jaune foncé.
Odour	Inodore.
Odour threshold	Data not available
pH	Solution aqueuse 5% = pH 8.6 - 9.8.
Melting point / Freezing point	Data not available
Initial boiling point	Data not available
Boiling range	Data not available
Flash point	Data not available
Evaporation rate	Data not available
Flammability	No
Lower flammable / Explosive limit	Data not available
Upper flammable / Explosive limit	Data not available
Vapour pressure	Data not available
Solubility	Très soluble dans l'eau. Insoluble dans l'alcool.
Vapour density	Data not available
Relative density	Data not available
Partition coefficient water/n-octanol	Data not available
Auto-ignition temperature	Data not available
Decomposition temperature	Data not available
Viscosity	Data not available

SECTION 10 - STABILITY AND REACTIVITY

Reactivity	Oxidizer: risk of fire in case of contact with combustible / organic substance. May react violently with incompatible substances.
Chemical stability	Stable under recommended storage conditions. If heated strongly, it emits toxic fumes.
Possibility of hazardous reactions	Oxidizer. Risk of fire or explosion if heated or crushed in presence of combustible or organic products.
Conditions of instability (Including sensitivity to shock / static discharge / vibration)	Avoid contact with incompatible materials.
Incompatible material	When pure, the product reacts with the following products: Reducing agents, Organic materials, Metallic powders. Hydrazine.
Hazardous decomposition products	Hazardous decomposition products formed under fire conditions. - Potassium oxides, Chromium oxides.

SECTION 11 - TOXICOLOGICAL INFORMATION

POTASSIUM CHROMATE

Routes of exposure	Ingestion, inhalation, skin and eyes.
Acute exposition effects / symptoms:	By exposure route below.
- Eyes	Severe irritation and burns that may cause permanent eye damage.
- Skin	Irritation and dermatitis. Exposed areas of the skin, including the fingers, the back of the hand and the forearms may present small ulcerations ("chrome sori").
- Inhalation	May be harmful if inhaled. Spasms, irritation and inflammation of the nose, throat and lungs. Edema of the larynx and bronchi. Chemical pneumonitis and pulmonary edema that can lead to death.
Acute toxicity (Ingestion)	Toxic if swallowed. Severe irritation of the digestive tract. Dysphagia, abdominal pain, gastrointestinal bleeding, cramps, diarrhea, melena, intense thirst, nausea and vomiting, tubular nephritis and possible hepatotoxicity, convulsions, circulatory collapse, unconsciousness, coma and can lead to death.
Chronic exposure effects / symptoms	Burning sensation, dermatitis, conjunctivitis, respiratory allergies, cutaneous reactions (chromium sori), kidney and lung damage, chest pain, cough, dyspnea, laryngitis, pharyngitis, headache, dizziness, ulceration and perforation of the septum nasal, metallic taste in the mouth, fatigue, digestive disorders (inflammation and ulceration of the gastrointestinal tract), weight loss and loss of appetite, seizures, nausea and vomiting. Liver damage. Skin sensitizer. Prolonged exposure to this product is likely to cause cancer. Is recognized as a carcinogen (class 1) by IARC.
DL50 (specify species and route of entry)	LD50 Oral - Rat - 7 - 168 mg/kg LD50 Oral - Mouse - 180 mg/kg LD50 Dermal - Rabbit - 2000 mg/kg.
CL50 (specify species and route of entry)	LC50 Inhalation - Rat - 4h - 83 - 99 mg/m ³ air

SUMMARY

Acute exposure effects / Symptoms:	By exposure routes below.
Ingestion	To our knowledge, the product has not been fully evaluated
Inhalation	To our knowledge, the product has not been fully evaluated
Skin	To our knowledge, the product has not been fully evaluated
Eyes	To our knowledge, the product has not been fully evaluated
Chronic exposure effects / Symptoms:	To our knowledge, the product has not been fully evaluated
ETA Mix (Estimated Acute Toxicity)	LD50 Oral: 737 mg/kg - Rat LD50 Dermal: > 5000 mg/kg - Rabbit LC50 Inhalation: 913mg/m ³ - 4h - Rat

SECTION 12 - ECOLOGICAL INFORMATION

Ecotoxicity	Potassium chromate: Toxicity to fish: LC50 - Pimephales promelas (fathead minnow) - 40 mg/L - 96 h Toxicity to daphnia and other aquatic invertebrates: EC50 - Daphnia magna (Water flea) - 15 mg/l - 48h Toxicity to algae: CE50- Nitzschia sp. - 0.26 mg/l - 72 h
Persistence and degradability	Data not available.
Bioaccumulative potential	Data not available.
Mobility in soil	Data not available.
Other adverse effects	An environmental hazard cannot be excluded in the event of unprofessional handling or disposal. Very toxic to aquatic life with long lasting effects.

SECTION 13 - DISPOSAL CONSIDERATIONS

Waste Disposal Method	Dispose of contents / container in accordance with local / regional / national / international regulations / or contact a specialist waste disposal company.
Contaminated Packaging	Dispose of as unused product.

SECTION 14 - TRANSPORT INFORMATION

UN Number	N/R
UN Proper shipping name	
Transport hazard class(es)	
Packing group	
Limited quantity index	
ERAP Index	
Special precautions	

SECTION 15 - REGULATORY INFORMATION

WHIMS CANADA	Skin corrosion/irritation - Skin irritation category 2 Serious eye damage/eye irritation - Serious eye damage category 1 Respiratory or skin sensitization - Respiratory sensitizer category 1 Acute toxicity - Oral category 4 Germ cell mutagenicity category 1 Carcinogenicity category 1 Respiratory or skin sensitization - Skin sensitize category 1
--------------	---

SECTION 16 - OTHER INFORMATION

Further information

The information in this document is based on the present state of our knowledge and is applicable to the product with regard to appropriate safety precautions. The above information is believed to be correct but does not purport to be all inclusive and shall be used only as a guide. It does not represent any guarantee of the properties of the product. Laboratoire MAT Inc. shall not be held liable for any damage resulting from handling or from contact with the above product.

Last Update: 4/14/2023