

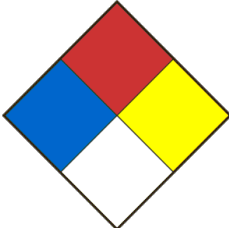


SAFETY DATA SHEET

SECTION 01 - PRODUCT AND COMPANY IDENTIFICATION

Product Identifier STEARIC ACID		Product Use Laboratory use	
Chemical formula C ₁₈ H ₃₆ O ₂		Product code SP-0155	Molar weight 284,48
Chemical name / Commercial name / Synonymous STEARIC ACID, ACIDE OCTADÉCANOÏQUE, OCTADECANOIC ACID, 1-HEPTADECANECARBOXYLIC ACID, PEARL STEARIC, CETYLACETIC ACID, STEAROPHANIC ACID, EMERSOL 132			
Supplier's name Laboratoire MAT		Address-Street 610, Adanac Street	
City Québec		Province Québec	
Postal code G1C 7B7	Internet www.labmat.com	Phone number 418-660-8666 Lun-Ven 8h-16h	
Emergency phone	418-660-8666 Lun-Ven 8h-16h		CENTRE ANTI-POISON DU QUÉBEC 800-463-5060
Date SDS 8/2/2024	SDS Prepared by Laboratoire MAT	E-Mail labmat@labmat.com	

SECTION 02 - HAZARDS IDENTIFICATION

Classification WHIMS / GHS	Not a hazardous substance according to WHMIS 2015
Other dangers	NFPA (Risk: 0=No risk; 1=Slight; 2=Moderate; 3=Signifiant; 4=Extreme)
	Health 1 Fire 1 Reactivity 0 Special danger

SECTION 03 - COMPOSITION/INFORMATION ON INGREDIENTS

Ingrédients (Dénomination chimique / synonymes)	Numéro CAS et tout identificateur unique	Concentration (%P/P)
Acide stéarique	57-11-4	<=100

SECTION 04 - FIRST AID MEASURES

Eye contact	Wash eyes with large amounts of water for at least 15 minutes while holding eyelids apart to rinse eyes. If irritation persists, seek medical attention.
Skin contact	Wash skin with plenty of water for at least 15 minutes. Remove soiled clothing. If irritation persists, seek medical attention.
Inhalation	Move the unwell person to the fresh air. If breathing is difficult, give oxygen. Consult a physician.
Ingestion	If the person is conscious, rinse the mouth with water. Never give anything by mouth to an unconscious person. Consult a physician.
Most important symptoms and effects (acute and delayed)	Main symptoms of high exposure: Irritation of the nose and throat. Eyes irritation. Ref. section 11.
Immediate medical attention and special treatment, if necessary	Treat according to symptoms. Show this sheet to the attending physician.

SECTION 05 - FIREFIGHTING MEASURES

Flammability	No
Suitable extinguishing media	Use water spray, alcohol-resistant foam, dry chemical or carbon dioxide.
Unsuitable extinguishing media	Data not available.
Hazardous combustion products	Hazardous combustion products formed under fire conditions: Carbon oxides.
Specific hazards of the dangerous product	Can form clouds of combustible dust. May react violently with incompatible products (Ref Section 10).
Special protective equipment and precautions for firefighters	Discard incompatible substances if this can be done without risk. Firefighters should be equipped with standard protective equipment, fireproof clothing, face mask, gloves, protective boots and, where appropriate, self-contained breathing apparatus.

SECTION 06 - ACCIDENTAL RELEASE MEASURES

Methods and materials for containment and cleaning up / Personal precautions, protective equipment	Evacuate personnel to safe areas. Ensure adequate ventilation. Use personal protective equipment. Pick up with a shovel or broom, taking care not to scatter dust. Dispose of residues in a container provided for the disposal of hazardous materials. Do not let product enter drains.
---	--

SECTION 07 - HANDLING AND STORAGE

Conditions for safe storage	Store in a cool, dry place. Keep container tightly closed and store away from heat, moisture, and incompatible products.
Methods of handling	Avoid formation of dust and aerosols. Provide appropriate exhaust ventilation at places where dust is formed. Avoid contact with the skin, eyes and clothes. Apply the usual standard hygiene rules: Wash your hands after use. Do not eat or drink during use.

SECTION 08 - EXPOSURE CONTROLS/PERSONAL PROTECTION

Workplace control parameters

Components	CAS-No.	Control	Value	Basis
STEARIC ACID	57-11-4	TLV, TWA, STEL	No data available	Canada. Alberta, Occupational Health and Safety Code (table 2: OEL)
		TLV, TWA, STEL	No data available	Canada. British Columbia OEL
		TLV, TWA, STEL	No data available	Québec. Regulation respecting occupational health and safety, Schedule 1, Part 1: Permissible exposure values for airborne contaminants

Ventilation	Use fan.
Respiratory	If permitted levels are exceeded, use NIOSH cartridge respiratory protection, or an air-supplied respirator.
Gloves	Handle with gloves. Suggested material: Nitrile. Butyle. Neoprene. The type, thickness and length of the glove must be chosen according to the use, the concentration of the product, as well as the duration of use. Replace gloves regularly for better protection.
Eyes	Safety goggles with safety shutters.
Shoes	Safety shoes.
Clothing	Labcoat. The type of protective equipment must be selected according to the concentration and amount of the dangerous substance at the specific workplace.
Engineering control	Have safety showers and eyewash stations in the workplace in case of an emergency and a ventilation system to maintain the level of concentrations in the air below the exposure limit values.

SECTION 09 - PHYSICAL AND CHEMICAL PROPERTIES

Physical state	Solid.
Appearance	Blanche.
Odour	Légère odeur évoquant le suif..
Odour threshold	20.00ppm
pH	Données non disponibles.
Melting point / Freezing point	67-72°C
Initial boiling point	361°C
Boiling range	Data not available
Flash point	113°C
Evaporation rate	Data not available
Flammability	No
Lower flammable / Explosive limit	Data not available
Upper flammable / Explosive limit	Data not available
Vapour pressure	1.0 mmHg @ 173.7°C.
Vapour density	9.8 (Air = 1.0)-
Relative density	0.845g/cm ³
Solubility	Insoluble dans l'eau. Soluble dans l'alcool, l'acétone, le chloroforme et l'éther.
Partition coefficient water/n-octanol	log Pow :8.23-
Auto-ignition temperature	395°C
Decomposition temperature	Data not available
Viscosity	12 mm ² /s à 70 °C.

SECTION 10 - STABILITY AND REACTIVITY

Reactivity	Acid product, reacts strongly with strong bases. May react violently with incompatible substances.
Chemical stability	Stable under recommended storage conditions.
Possibility of hazardous reactions	Stable under normal conditions. May react violently with incompatible substances.
Conditions to avoid, including static discharge, shock or vibration	Avoid excessive heat.
Incompatible materials	Strong oxidizing agents (nitric acid, perchloric acid, peroxides, chlorates and perchlorates), reducing agents (potassium, sodium, metal hydrides), bases, heat and moisture.
Hazardous decomposition products	Hazardous decomposition products formed under fire conditions. Carbon oxides.

SECTION 11 - TOXICOLOGICAL INFORMATION

STEARIC ACID

Routes of exposure	Ingestion, inhalation, skin and eyes.
Acute exposition effects / symptoms:	By exposure route below.
- Eyes	Irritation.
- Skin	Irritation.
- Inhalation	To our knowledge, the product has not been fully studied.
Acute toxicity (Ingestion)	Irritation of the mucous membranes. Gastrointestinal disorders, cramps, diarrhea, headache, dizziness, nausea and vomiting.
Chronic exposure effects / symptoms	Burning sensation. Cough, gastrointestinal discomfort, irritating effects.
DL50 (specify species and route of entry)	LD50 Oral - Rat - 2000 mg/kg LD50 Dermal - Rabbit - >5000 mg/kg
CL50 (specify species and route of entry)	LC50 - Inhalation - Data not available.

SECTION 12 - ECOLOGICAL INFORMATION

Ecotoxicity	LC50 - Lepomis macrochirus - Static test: 4500 ug/L- 24 h. LC50 - Oncorhynchus mykiss (rainbow trout) - Flow-through test - 3350 ug/L - 24 h.
Persistence and degradability	Data not available.
Bioaccumulative potential	Data not available.
Mobility in soil	Data not available.
Other adverse effects	May be harmful in the long term for aquatic organisms.

SECTION 13 - DISPOSAL CONSIDERATIONS

Waste Disposal Method	Dispose of contents and container in accordance with local, regional and national regulations, or contact a specialist waste disposal company.
Contaminated Packaging	Dispose of as unused product.

SECTION 14 - TRANSPORT INFORMATION

UN Number	NR
UN Proper shipping name	
Transport hazard class(es)	
Packing group	
Limited quantity index	
ERAP Index	
Special precautions	

SECTION 15 - REGULATORY INFORMATION

WHIMS CANADA	Not a hazardous substance according to WHMIS 2015
---------------------	---

SECTION 16 - OTHER INFORMATION

Further information

The information in this document is based on the present state of our knowledge and is applicable to the product with regard to appropriate safety precautions. The above information is believed to be correct but does not purport to be all inclusive and shall be used only as a guide. It does not represent any guarantee of the properties of the product. Laboratoire MAT Inc. shall not be held liable for any damage resulting from handling or from contact with the above product.

Last Update: 8/2/2024